


Centrifuge CM-7S

Advanced technology



- SkySpin™ patented motor controller.
- Sensorless Accu-Logic architecture
- Digital speed and balance controller.
- Increased durability.
- Smooth, noiseless and non-heating operation.

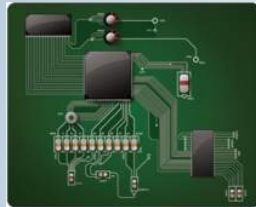

Advanced interface



Major control buttons: "Start/Stop" and "Open" are comfortably located right underneath the lid handle. No matter right or left handed user can always reach the control buttons with a thumb.



Large colorful LCD display makes navigation and settings very simple and fast. All the settings are displayed simultaneously which is very convenient.



Easily accessible settings menu allows additional fine tuning. User can adjust display brightness and sound volume to insure maximum comfort.



Low noise and heating



Lid lock



Service free drive



Rotor imbalance sensor

TECHNICAL SPECIFICATIONS

Rotor rotation speed	RPM	100 - 3500
Timer	min	1 - 99
Maximum centrifugal force	RCF	2300
Brake levels		6
Working range of temperatures	°C	from +10 up to +40
Relative air moisture at 20 °C	%	80
Voltage and frequency	V,Hz	100 - 240V; 50/60Hz
Maximal consumable power	W	250
Size of the device (L x W x H)	mm	430 x 410 x 220
Weight	kg	12.6